



H360

3 Part Diff Hematology Analyzer





22 Parameters Enhanced patient care



Supports geriatric and paediatric samples



Convenient patient data review with customizable reports

Making Automation Affordable for labs Everywhere

H 360

3 Part Diff Hemetology Analyzer

Erba H360 is the new generation 3 Part Diff Analyzer equipped with advanced clinical parameters for efficient clinical diagnosis.


KEY FEATURES



RBC / Platelet electrical impedance

22 Parameters including **PLCC & PLCR**

• Electrical impedance measures change in electrical resistance produced by a blood cell suspended in conductive diluent as it passes through an aperture

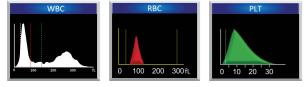

- P-LCR, P-LCC, PDW-SD & PDW-CV aids in differential diagnosis of cases with abnormal platelet counts
- RDW-SD, RDW-CV & MCV aids in screening for causes of anaemia




Convenient

- 10.4" large coloured LCD touchscreen with graphic icons
- Single page view of full report
- ADDM Unique Auto Diluent Dispensing **M**echanism
- RIMS Real time Inventory Management System with RFID






- Excellent short and long-term precision within 2%
- Highly accurate results correlating with microscopic examination & clinical findings
- Specific flagging messages to alarm for microscopic slide review

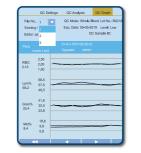



- Just 9 µL sample volume: All age groups
- Low blood sample volume

- Consolidated patient reports can be generated for customized parameters such as blood type, Rh blood group, ESR, reticulocyte, CRP along with other reportable parameters
- QR code scanner facilitates direct inputs of QC values

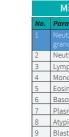

- Lab mean
- Traceability of controls and calibrator runs

Assured Quality with

NABL compliance

• Traceability of consumable lots



Advanced Software

• Provision to include microscopic investigations such as bands segment, blast, etc. in the printed report

icroscopic Exam. Settings			
imeter Name	Code System		
		*	
rophilic band granulocyte			
ohocyte			
ocyte		_	
nophil			
phil			
macyte			
ical Lymph		¥	
-			

Ordering Information

Product Code	Description	Pack size (ml)
166122	Erba H360 analyzer	-
121740	H 360 Dil	20,000
121864	H 360 Lyse	1 × 500
121741	H 360 Lyse	3 x 500
121938	H clean	1 x 50
121939	H clean	4 x 50

Erba Lachema s.r.o.

Karásek 2219/1d, 621 00 Brno, Czech Republic Phone: +420 517 077 111 ✓ sales@erbamannheim.com

www.erbamannheim.com



Transasia Bio-Medicals Ltd. An ISO 13485:2016 / 9001:2015 Certified Company Transasia House, 8, Chandivali Studio Road, Chandivali, Andheri (E), Mumbai - 400 072 T: (022) 4030 9000 F: (022) 2857 3030

✓ responses@transasia.co.in